


# F&P AIRVO 2

AIRVO™ 2 humidifier with integrated flow generator

The AIRVO 2 generates high flows of warmed and humidified respiratory gases, delivered to the patient through a variety of nasal, tracheostomy and mask interfaces.

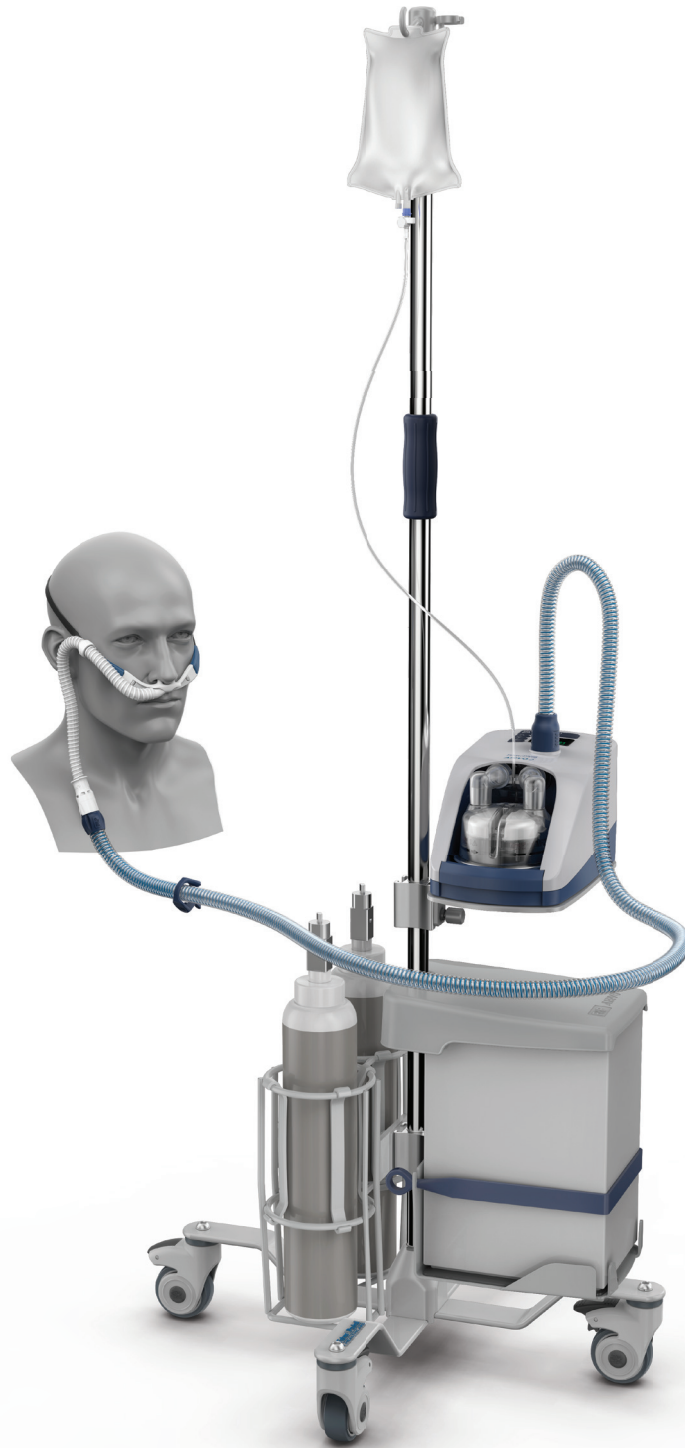

# AIRVO 2 | FEATURES

With its integrated flow generator and innovative oxygen delivery system, the AIRVO 2 provides an effective, versatile solution for the delivery of oxygen therapy.

## High Performance

- Designed to exceed the ISO standard for respiratory humidification
- Up to 60 L/min of supplemental oxygen


## Versatile

- Variety of Optiflow patient interfaces for upper and bypassed airways
- Wide flow range: 2-60 L/min

## Set up and use

- Graphical display menu for settings
- Integrated O<sub>2</sub> mixing
- Inbuilt O<sub>2</sub> sensor
- No probes or external air supply required

## ACCESSORIES




Optiflow™+ Nasal Cannula

- OPT942 (small, 20-pack)
- OPT944 (medium, 20-pack)
- OPT946 (large, 20-pack)




Optiflow™ Junior Nasal Cannula (20-pack)

- OPT316/OJR416 (infant)
- OPT318/OJR418 (pediatric)




Optiflow™+ Tracheostomy Interface

- OPT970 (20-pack)



Optiflow™+ Mask Interface Adapter

- OPT980 (20-pack)

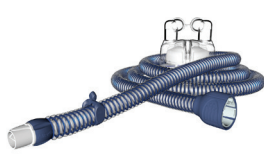

Tube & Chamber Kit

- 900PT561 (10-pack)




Disinfection Kit

- 900PT600



Cleaning Sponge Stick

- 900PT602 (20-pack)

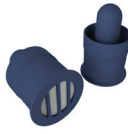

Disinfection Filter

- 900PT601 (2-pack)



Air Filter

- 900PT913 (2-pack)



Oxygen inlet kit

- 900PT422



Pole mounting tray

- 900PT405



Hospital Pole Stand

- 900PT421



UPS mounting kit

- 900PT411



Oxygen bottle holder

- 900PT427

AIRVO 2 HUMIDIFIER WITH INTEGRATED FLOW GENERATOR	
Intended use	The AIRVO 2 humidifier is for the treatment of spontaneously breathing patients who would benefit from receiving high flow warmed and humidified respiratory gases. This includes patients who have had upper airways bypassed. The flow may be from 2 - 60 L/min depending on the patient interface. The AIRVO 2 is for patients in hospitals and long-term care facilities.
Humidity output	Complies with ISO 8185:2007 standard
Flow settings	10 to 60 L/min (in increments of 5 L/min) (default) 2 to 25 L/min (in increments of 1 L/min) (Junior mode)
Patient interfaces	F&P Optiflow™+ nasal cannulas F&P Optiflow™+ tracheostomy interfaces Masks with standard 22 mm medical taper
Models	PT101AZ / PT101EE / PT101EW / PT101UK / PT101US / PT101AN / PT101JP / PT101KR
Dimensions	295 mm x 170 mm x 175 mm (11.6" x 6.7" x 6.9")
Weight	2.2 kg (4.8 lb) - unit only 3.4 kg (7.5 lb) - packaged in bag incl. accessories
Electrical ratings	50-60 Hz 100-115 V - 2.2 A (2.4 A max) 220-240 V - 1.8 A (2.0 A max)

Some products may not be available in your country. Please contact your local Fisher and Paykel Healthcare representative.